



User Manual

Mandala Breeze



Please read the manual before using the tools.

Version 1.0
Updated February 26, 2025

1. General Notes

1.1 Product Feature and Intend Use

Mandala Breeze is a sleek, portable wellness device that helps **restore balance and reduce stress** using magnetic fields, light therapy, and gentle vibrations. Its quantum and bifilar coils fine-tune your body's energy flow, promoting relaxation and well-being.

Pair it with the Mandala Breeze app, which allows you to upload and create custom presets for personalized energy balancing. The app also includes basic pre-set programs for immediate use. Additionally, by integrating **Mandala Scan**, you can analyze your body's energy state and identify the best coherence preset to upload for optimal results.

Your **Mandala Breeze** has been **quality-tested** and is ready to support your wellness journey anytime, anywhere.

1.2 Safety Warnings

- Check all accessories upon unboxing.
- Read and follow all safety instructions in the manual.

1.3 Warranty

- Your Mandala Breeze includes a one-year warranty from the date of delivery.
- If you experience any issues, contact customer support immediately.
- The warranty does not cover damage due to:
 - Misuse or neglect
 - Unauthorized modifications or disassembly
 - Customs/import fees or indirect damage
 - If repairs are needed due to misuse, all costs (including shipping) will be the user's responsibility.

1.4 Important Precautions

- Do not use the device if it appears damaged. Contact support.
- Ensure the correct voltage for your country before use.
- Charge for at least 1 and half hours before the first use.
- Do not disassemble or modify the device.
- Keep dry – avoid rain, high humidity, heat, and dust.

1.5 Maintenance

1. If unused for a while, charge the device for 2 hours before use.
2. Use soft cloth or water-based glass cleaner for maintenance.
3. Keep the device clean to ensure proper performance.

1.6 Product Includes

1. Mandala Breeze
2. USB Type-C cable
3. Package box




2. Technical Data

2.1 Specifications

Product name:	Mandala Breeze
Color:	White
Material:	ABS plastic
Input Voltage:	DC 5V / 2A Max
Charging Interface:	USB Type-C
Battery Type:	3.7v / 500 mAh Li-Po Battery
LED Output:	LED RGB 3528 / UV-C LED
Coil Output:	Bifilar Coil / Quantum Coil
frequency:	0 Hz – 2 MHz
Bluetooth:	BLE Bluetooth Standby Mode
Length:	45 mm.
Width:	45 mm.
Height:	18 mm.
Weight:	320 grams

2.2 Overview and Function

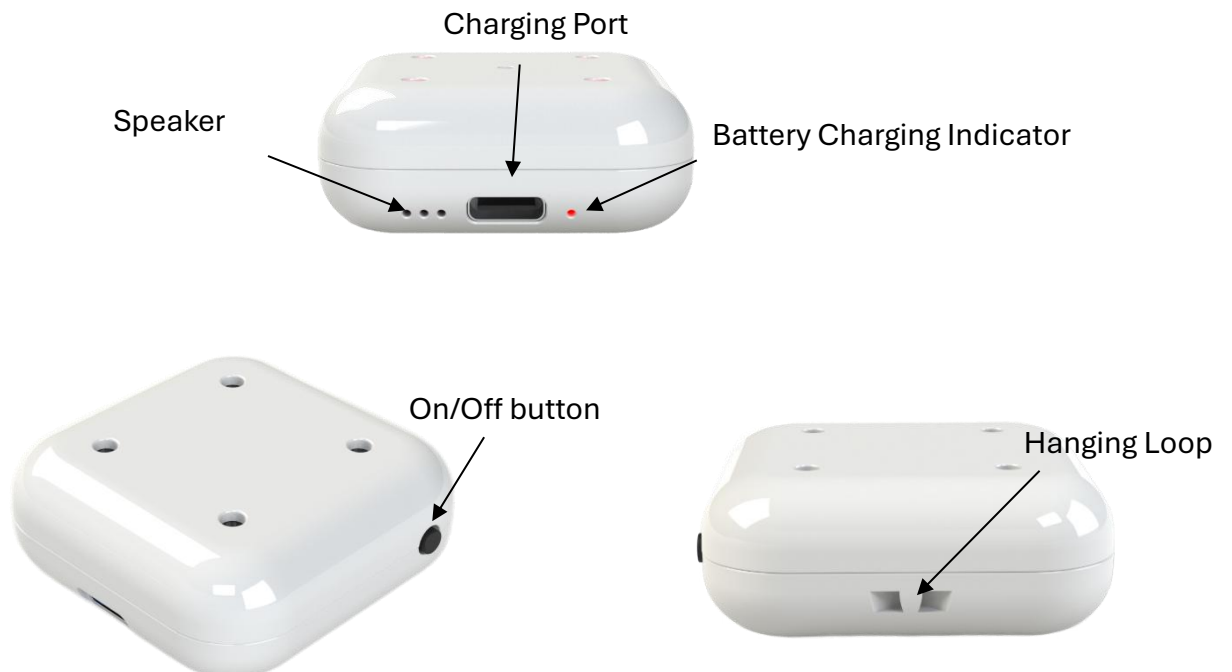

2.3 Charging Setup

1. To charge the Mandala Breeze, connect the USB-C cable to the device and a power adapter, then plug it into a power source.

2. Charging takes about 1 and half hours from an empty battery and lasts over 4 hours. *(Battery capacity can be reduced.)*

3. The red light will blink while the battery is charging and will stop once the battery is fully charged.




To protect the battery life span, please don't overcharge the battery.

3. Mandala Breeze Application


3.1 Application Setup

The Mandala Breeze app is available on iOS and Android. Download it from the App Store for iOS or Google Play Store for Android. For users in China, please download it from <https://light-mandalas.com/downloads>

App Store (iOS)

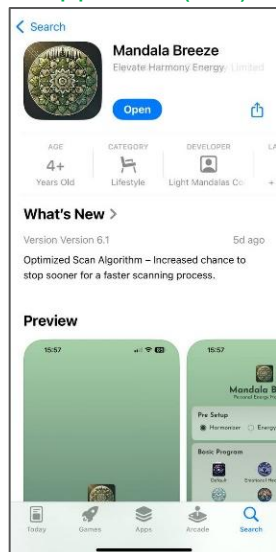


Play Store (Android)

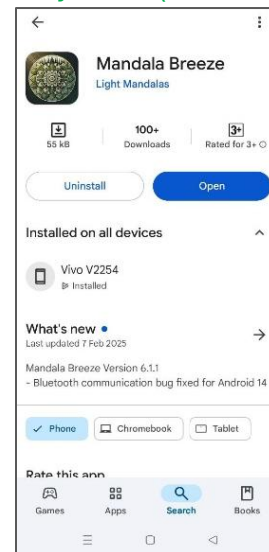

3.2 Pairing Bluetooth with Application

Tap the "Settings" icon at the bottom right, then select "Mandala Breeze Setup." The app will search for nearby devices. Choose the displayed serial number to connect via Bluetooth.

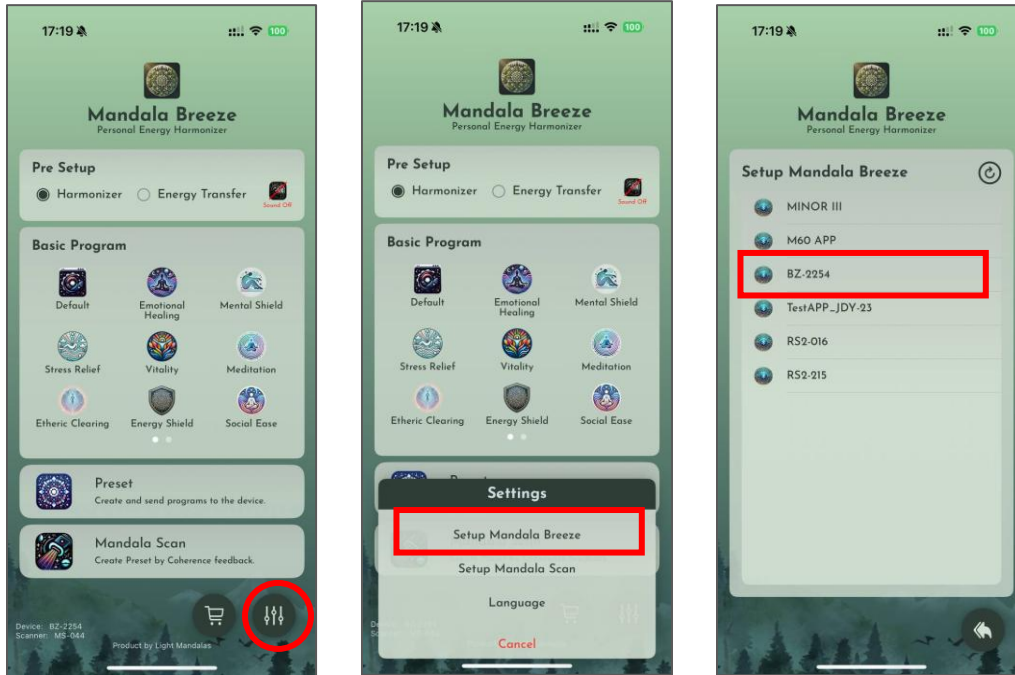

3.3 Pre-Setup

Harmonizing Mode:
Balances and unifies different energies, reducing conflict and creating a calmer, more stable atmosphere.





Sound On / Off

Torsion Mode:
Generates subtle, vortex-like fields that hold energies and encourage a fresh, flowing energy from the past to present.

3.4 Basic Program

After setting up the app, select your preferred program under "Basic Program" and tap "Upload" to upload it.

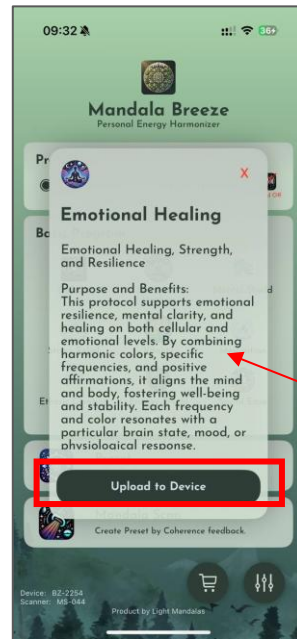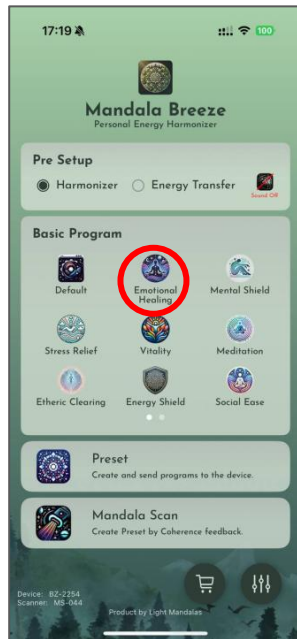



Description and Basic Instruction

3.5 Preset

"Preset" allows users to customize frequency settings using the library or their own frequencies. Users can combine frequencies to create a preset.


3.5.1 Frequency Library

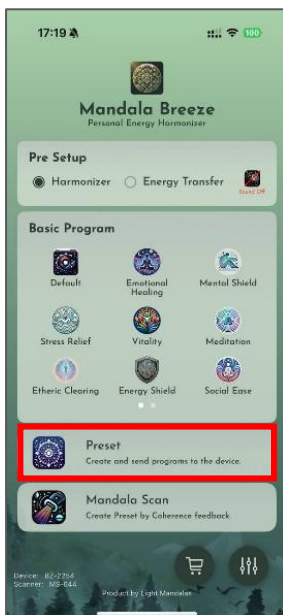

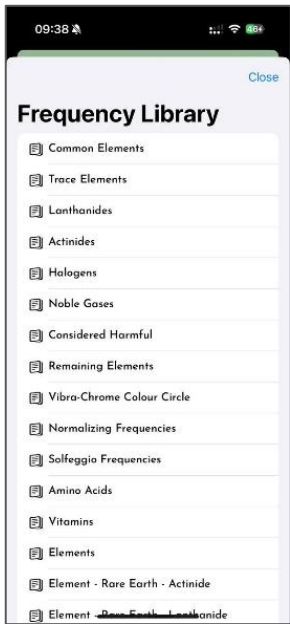


1. Give the name for the preset.

Select the frequency that you prefer to add to the preset, then click "Confirm".

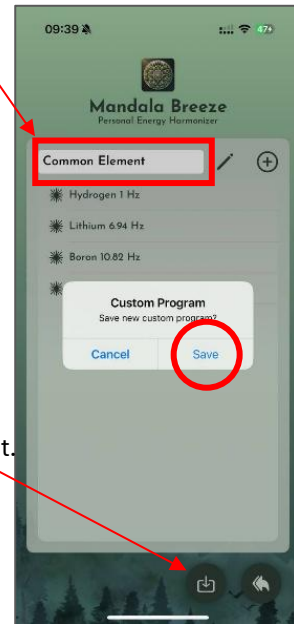


2. Click this icon to save the preset.

3.5.2 Insert Frequency

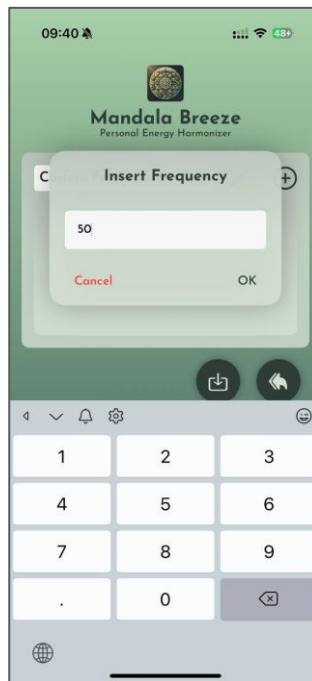


1. Give the name for the preset.

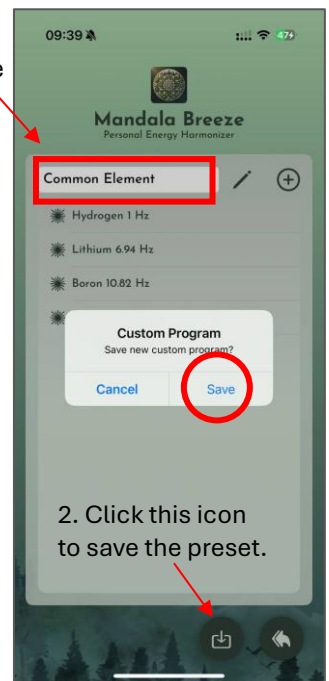


2. Click this icon to save the preset.

3.6 Mandala Scan

To use Mandala Scan, first pair it via Bluetooth. Tap the "Settings" icon, select "Mandala Scan Setup," and choose the displayed serial number to connect.

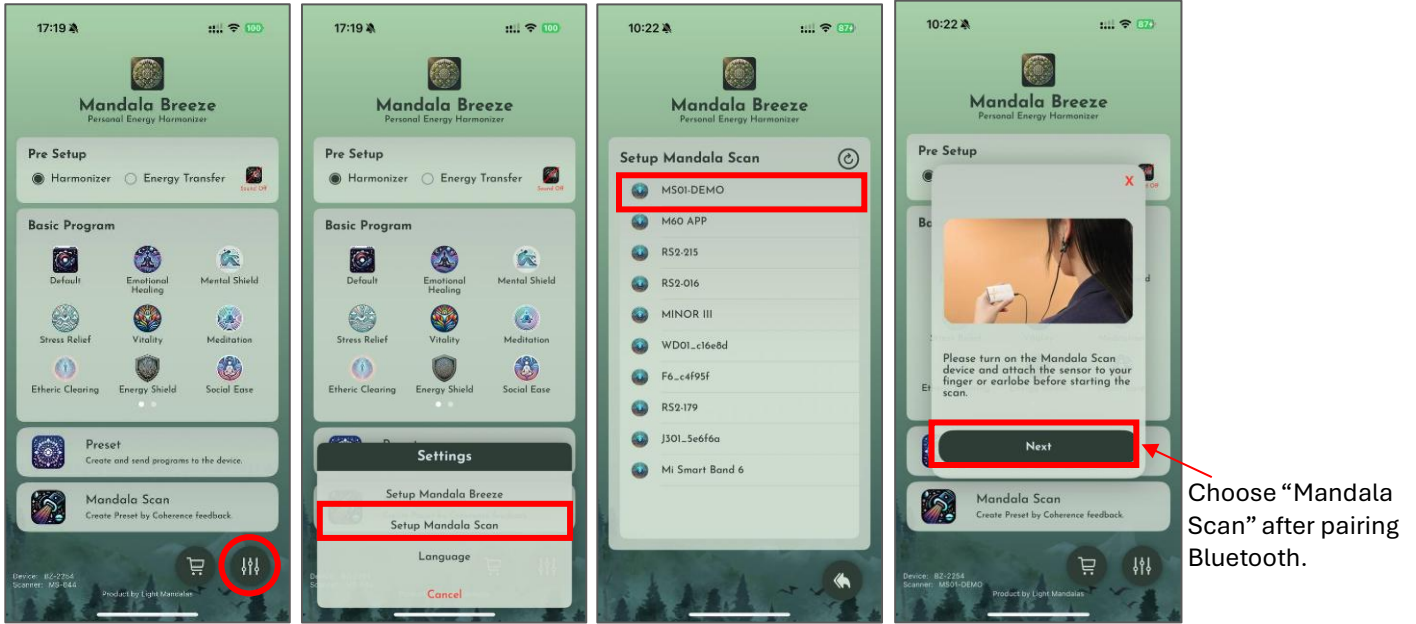


3.6.1 Smart Scan

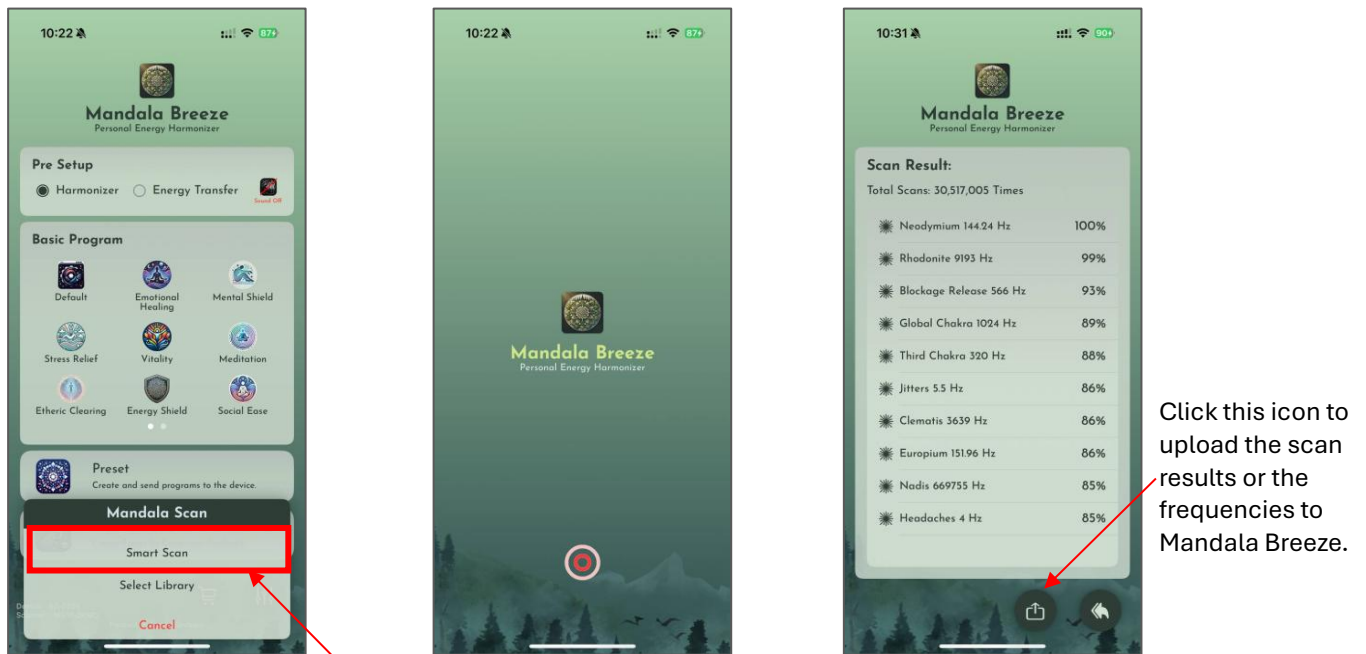



Smart Scan analyzes your body and identifies the frequencies that best support your needs at this moment.

Click this icon to upload the scan results or the frequencies to Mandala Breeze.


3.6.2 Select Library

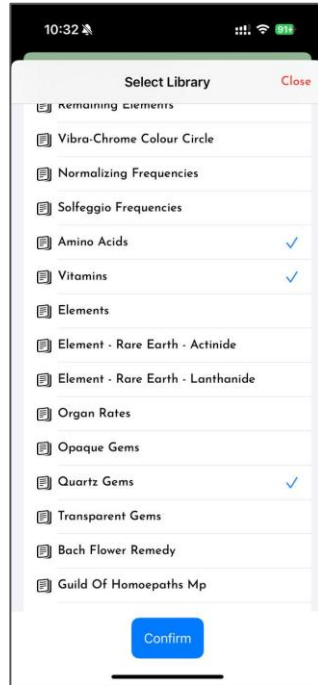

Select the library you want to scan, then click 'Confirm' to begin scanning.

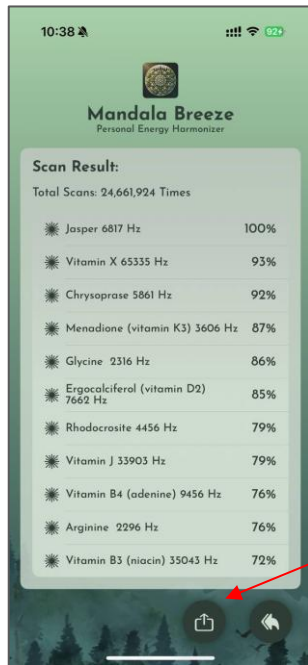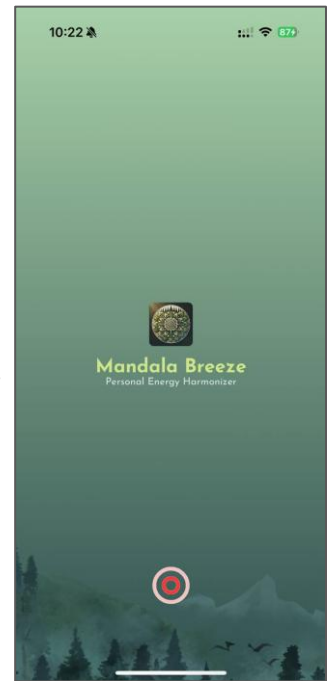

Click this icon to upload the scan results or the frequencies to Mandala Breeze.

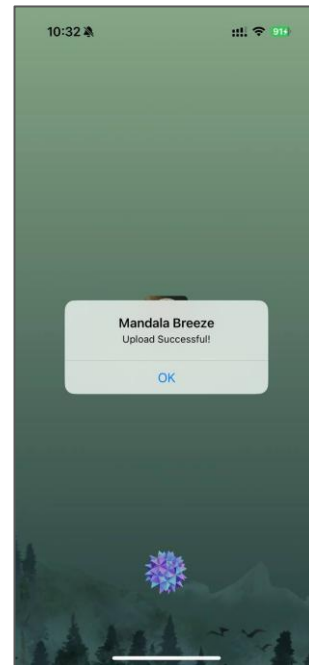

Failures Remedies

Device cannot turn on

- The battery is low. Please charge the device.

Disclaimer

Light Mandalas products are designed as energy support tools for relaxation, stress management, and well-being. They are not medical devices and do not replace professional medical treatment or advice.

Any experiences or testimonials shared regarding the use of our products reflect individual responses and are not guaranteed outcomes. Since personal energy and emotional states vary, users may experience different effects, including heightened sensitivity or discomfort. By using our products, you acknowledge that all outcomes are your sole responsibility.

Light Mandalas Co., Ltd. (Thailand) is not liable for any unintended effects, discomfort, or injuries that may arise from product use. If you have any medical concerns, please consult a healthcare professional before using our products.

For any questions, please reach out to us for further support.

Light Mandalas Co., Ltd. (Thailand Headquarters):

- Light Mandalas Website: www.light-mandalas.com
- Email: info@light-mandalas.com
- Instagram: [@mandalawave9](https://www.instagram.com/mandalawave9)
- YouTube Channel: Light Mandalas Official
- TikTok: Light Mandalas
- WeChat: Light Mandalas 光之曼陀罗
- Technical Support Email: cservice@light-mandalas.com